

MK-P29

Emergency Trolley



*Better care
Better comfort!*

MK-P29

EMERGENCY TROLLEY

DESCRIPTION

The durable stainless steel countertop and the automatic return function of the drawer make the mk-p29 a popular choice in emergency cart. The perfect configuration and smooth operation experience can greatly improve the work efficiency of doctors and become an indispensable product in emergency scenarios.

FEATURES

- 1.Made of aluminum alloy frame and ABS plastic structure, light and firm
- 2.Stainless steel countertop ,antibacterial and corrosion resistant
- 3.5 drawers with dividers and card slots,has automatic return function
- 4.Detachable storage box, garbage can,basket,convenient to use.
- 5.With IV pole, Defibrillator shelf ,CPR board and Oxygen tank holder
- 6.Center-lock system with key or disposable sealing
- 7.Four dia 100mm noiseless castors with 2 brakes
- 8.Ergonomically designed handle for easy movement

STANDARD SPECIFICATION

Overall size	750*360*960mm
Countertop size	42*57cm
Cpr board size	41*53cm
Extended working surface size	32*37cm
Structure	Aluminum alloy frame and abs plastic structure
Castor	Dia 100mm




Extended working surface

Telescopic shelf which can increase the working surface.




Locking system

Central locking system with key or disposable seal.

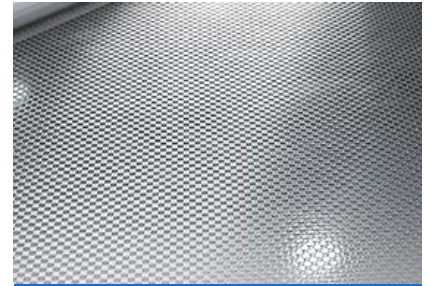

Countertop

ABS engineering plastic with stainless steel countertop, antibacterial and corrosion resistant.

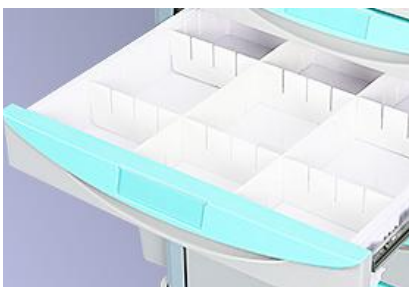

Drawer

5 drawers with dividers inner and card slot outside, also have automatic return function, which is convenient to use.



Storage boxes

Big storage space, convenient to use.



Garbage can

Garbage cans with different colors facilitates garbage sorting.



Basket with sharp boxes

Storage for sharp needles.



Power outlet

Power outlet for equipment using, normal length 3m.



CPR board

CPR board for patient rescue.




IV pole

IV pole with 4 hooks, height adjustable.



Defibrillator shelf

Angle and height adjustable, for equipment storage.

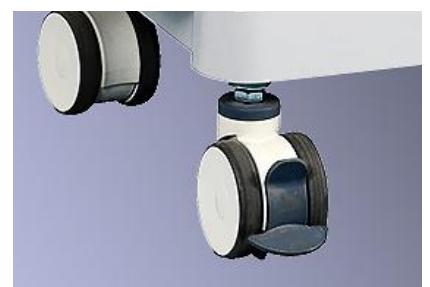

Castor

Four diameter 100mm noiseless castors with 2 brakes.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US