

YA-D3-5

Ultra Low Electric Hospital Bed



*Better care
Better comfort!*

YA-D3-5

ULTRA LOW ELECTRIC HOSPITAL BED

DESCRIPTION


YA-D3-5 ultra-low hospital bed has an Auto Contor function that simultaneously raises both the head and the knees, residents stay comfortably in place while being better positioned for reading, watching TV, or talking to family members. The simultaneous movement also helps prevent painful skin shearing that might otherwise occur during repositioning. It is the ideal bed with its stability, easy controls, and smooth raising and lowering of the head, feet, and overall height of the bed.

FEATURES

- Ultra-low with the minimum height at 260mm
- Cold rolled steel frame with powder coating treatment
- Electric motorize bed with smooth and control operation by 3 actuators
- Remote hand control
- Auto contour backrest and kneerest adjust simultaneously
- Aluminum alloy material foldable type side railing with 6 bars
- Detachable head/foot board
- Electric CPR function
- With 4 corners bumpers
- 4 provisions for IV pole
- 2 drainage bag holders


STANDARD SPECIFICATION

Overall size	L2050*W920*H260mm
Height adjustment range	260-560mm
Cators	125mm
Back board	0-75°
Knee board	0-40°
Loading weight	200kg

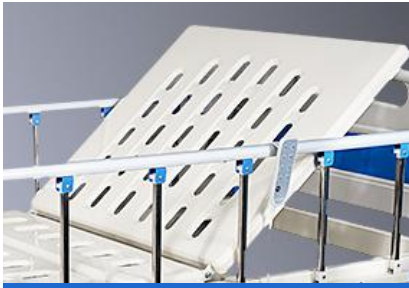

Back board

Back board can be adjust from 0 to 75 degree.

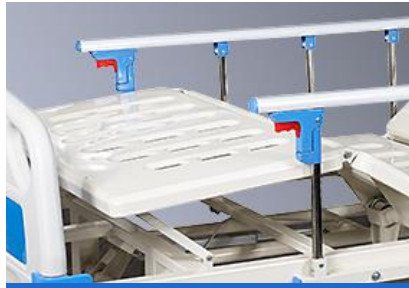

Knee board

Knee board can be adjust from 0-40 degree.



Side railing

Aluminum alloy fold-able type side railing with 6 bars.



Wheels and locking

Dia 125mm double side mute wheels, swivel type with independent locking.



Bump cover

There are 4 bump cover in each corner, it will protect bed well.



Detachable head/foot board

Head/foot board can be removable easily and shortly due to the detachable design.



Remote hand control with electric CRP

Nurse can adjust the bed via the remote hand control for different position.



Platform

Stamping platform with ventilator hole for skin breath.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US