

YA-D5-11

Electric ICU bed



*Better care
Better comfort!*

YA-D5-11

ELECTRIC ICU BED

DESCRIPTION

Recovering patients require a nurturing environment where they can recuperate with peace of mind. Placing particular emphasis on safety, ya-d5-11 electric bed is designed to mitigate accidents, In addition to safety, the beds are also designed for usability, easy cleaning, and easy maintenance.

FEATURES

- Four split side rail, drop down.
- Hand remote control.
- Detachable head/footboard with easy lock.
- Central brake system with a wheel diameter of 125mm.
- Bumper caster on four corners.
- Auto-counter to provide the patient from sliding onto foot end.
- Abs articulated soft joint for the platform.
- With angle indicator embedded in the railing for back and trendelenburg adjustment (optional).
- Manual cpr on both sides.
- Incorporated railing control (optional).
- Weighing scale with led display (optional).

SPECIFICATION

Overall size	2200*990*350-700mm
Recommended mattress size	1950*900mm
Height adjust	350-700mm
Back adjust	0° to 75°
Knee adjust	0° to 40°
Trendelenburg	0° to 12°
Reverse trendelenburg	0° to 12°
Load capacity	225kg



Side rail control

Incorporated railing control from both inside and outside of the railing panel, vastly convenient for both patient and doctor, and easy for hospital management.



Remote control

The hand control with clear icon is beneficial for medical staff which can control the bed position from a remote distance.



Emergency shut down button

Everything will stop once press the emergency shut down button if any mistake operation.



CPR function

Release the back-rest for emergency CPR.



Central locking

The bed will be locked and unlocked by the foot pedal, it is central locking system.




Detachable head&foot board

Both the head and foot board is detachable, and tighten the knob to secure the headboard and footboard.




Angle indicator

Easy-to-read angle indicators located on the.

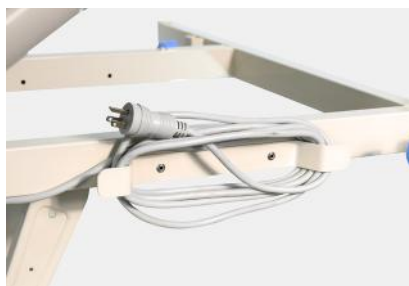

Wrap the power cord

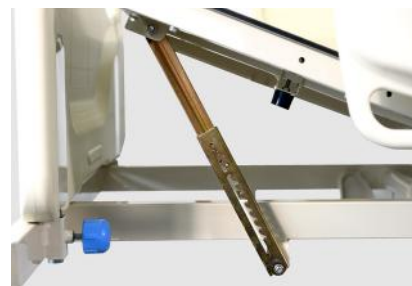

Foot rechat

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US