

# YA-D5-2

Electric hospital ICU bed



*Better care  
Better comfort!*

# YA-D5-2

## / ELECTRIC HOSPITAL ICU BED

YA-D5-2 defines a new standard for a popular ICU bed, the full length and streamline side rail, fashionable and ergonomic design, reinforced X shape structure with strong stability. It is an ideal choice for the major ICU room setup.

### Main Features

- Four section HDPE platform
- Full length guardrail with safety gap
- Control panel on both inside and outside of the guardrail
- Thoughtful handle design on the inside of the guardrail
- Removable ABS head and foot board with lock device
- Electric and manual CPR
- Auto-Counter
- One button cardiac chair position
- Locking function for control panel
- Battery backup (optional)
- Weighing scale (optional)
- Emergency button (optional)
- Night light
- Four provision for IV pole

### Technical parameter

Overall size	2240*1060*490-790mm
Lying surface size	1980*900mm
Back section	0-70°±5°
Leg section	0-40°±5°
Height adjustment	490-790±10 mm
Trendelenburg/Revere- Trendelenburg	0-12°±3°
Safe working load	250kg
Castor	Dia 125mm

## Main details



**Height adjustment:** 490-790±10mm



**Trendelenburg:** 0-12°±3°



**Reverse trendelenburg:** 0-12°±3°



**Back section:** 0-70°±5°

**Leg section:** 0-40°±5°

Back and foot section can be adjusted individually , also has auto-counter function



**Cardiac Chair Position**



### Side rail

Full-length side rails with safety gap, comply with IEC 60601-2-52 standard, tuck-away type, adjusted through the buckle.



### MS perforated platform

MS perforated platform with ventilation hole, molded.



### Nurse control

Nurse control on the foot end panel, specially has one button cardiac chair position.



### Hand remote controller




### Guardrail control system



### Control panel and handle

Control panel and handle on the inside of the guardrail are convenient for patient operation.

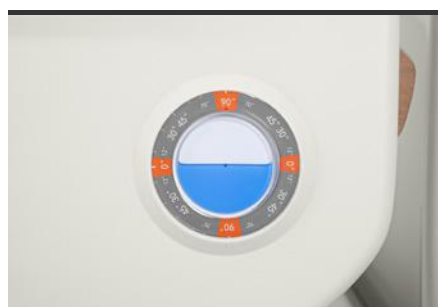


### Angle display

Angle display on the side rail, can indicate the adjustment angle of the bed.



### Manual CPR function

Manual CPR function on both sides of the bed for quick release.

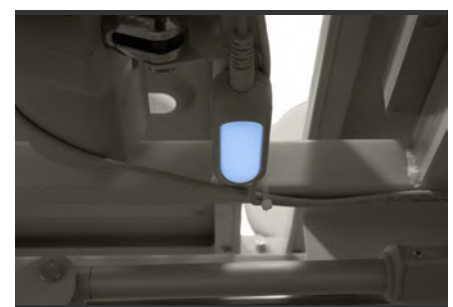


### Night light

The night light is soft and can provide illumination for patients and medical staff at night.



**Bumper castors**

Four bumper castors to protect the bed during transport.

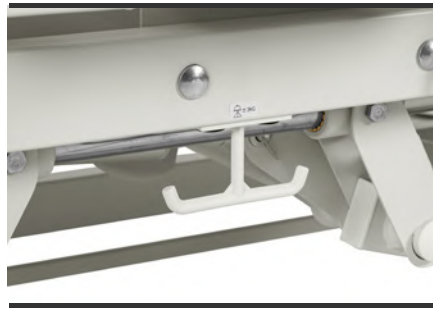

**Urine bag hook**

Urine bag hook under both sides of the bed.



**Mattress holder**



**IV pole**



**Bedside card on foot panel**



**Emergency stop button**

Removable head and foot board with lock device.



**Castors**

Four dia 125mm castors , central brake pedal on foot end of the bed.

## ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

**Pls scan the QR code for more info**

**info@medik-medical.com | +86 139 6223 1292**

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US