

YA-D5-8

Electric Hospital Bed



*Better care
Better comfort!*

YA-D5-8

ELECTRIC HOSPITAL BED

DESCRIPTION

The YA-D5-8 Five Function Electric Hospital Bed has extensive adjustability to make it as easy as possible for patients to receive high-quality care. Hand remote control operation is smooth. side rails work well to protect patients from falls. It is the best choice for hospital wards.

FEATURES

- Bed Frame Electrophoresis and powder coating
- Bed platform MS perforated platform with ventilation hole, molded
- Hand remote controller of all electrically operated functions and movements
- Handrails Collapsible aluminum railing with 6 bars
- Head and Footboard ABS material, detachable
- Castors: Central brake system with wheel diameter of 125mm
- With a 5th wheel for direction

TECHNICAL PARAMETER

Model	YA-D5-8
Overall size	2100*960*550mm (82.5" x 37.7" x 21.5")
Height adjustment	550~750mm (21.5" – 29.5")
Back board	0~75°
Knee board	0~40°
Trendelenburg	0~12°
Reverse trendelenburg	0~12°
Castor	125mm (5")
Loading capacity	200kg(440 lbs)

PRODUCT DETAILS

Hand remote control

The clear icon of the hand control is beneficial to the operation of staff.



Side rail

6 bar Aluminium alloy side rails, pushing button operates fold and open.

Option: PP side rails



5th wheel

The bed has 5th direction wheel provides the stretcher with ideal manipulation and stability.




Castor

125mm central control castor brake located at the foot of the bed.



Bumper castor

The bumper castors are on the four corners, which can reduce the impact of accidentally bumping into the wall during the movement.

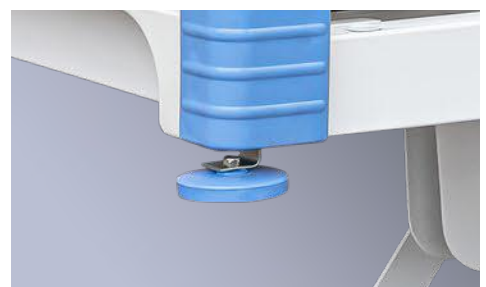

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US