



YA-DH3-1

Home care bed

Home care bed

YA-DH3-1

YA-DH3-1 is an electric bed with scissor action type design for home care and nursing care center. Full length protection side railing minimizes the risk of patient falling down, while multi functionality of the bed also meets various needs of patient comfort as well as convenience for medical care staff. Outside wooden style frame gives an home like environment for patient to keep them in good psychic condition.



Main Features

- Electrically operated profiling bed for community and nursing home use.
- Safe work load (SWL) 180Kg.
- Wooden head and footboard with aluminum columns around.
- Adjustable self-locking side rails providing additional safety which can be released at both head and foot ends.
- Remote handset control.
- High-quality castors (castor diameter: 75mm).
- Electric CPR.

Main details

YA-DH3-1

Height adjustment

The overall stroke of the home care bed is 360mm. When it is at the lowest position of 360mm, it is an ideal sleeping position, and such height will help use to exit bed easily. When it reaches the highest position of 760mm, it is convenient for medical examination.



Germany Dewert control system

Total 3pcs linear heavy duty type actuator will guarantee different position adjustment and safety loading 180kg



Bed platform

Sturdy metal welded platform surface with mattress holder, 4 section board.


Main details

YA-DH3-1



Side railing control button

The full-length side rail prevents the patient falling down, and the height of the guardrail can be adjusted by pressing button.

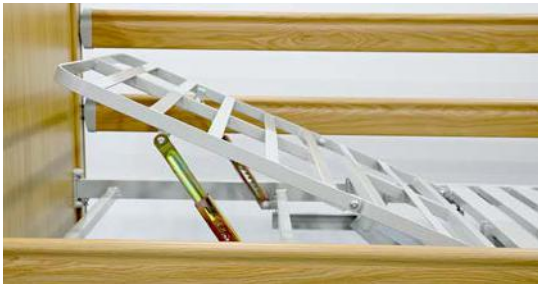

Calf board

Calf board can be adjusted by racks separately, which can raise level the feet, promote blood circulation if necessary.



Remote control with electric CPR

Hand Remote control with electric CPR function, millions of secondary button service life, clear function identification, easy operation.



Castor

Unique eight $\Phi 75\text{mm}$ swivel casters, high support, impact resistance, not easy to break.

Technical parameter

YA-DH3-1

Model	YA-DH3-1
Overall size	2120*1030*360mm
Lying size	1900*900 mm
Backrest adjustable	0-75°
Footrest adjustable	0-40°
Height adjustable	360-760mm
Castor size	75mm
Loading capacity	180kg

Options

- Lifting pole with triangle
- Urine bag stand
- Battery backup
- Extension kit
- Longer railing

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US