



YA-EC-U04

**Electric Cardiac Ultrasound
Examination Table**

YA-EC-U04

Unique cardiac ultrasonic scanning position, automatic paper rolling and auto-correction system make the YA-EC-U04 as an outstanding exam table specially for complicated diagnosis procedure. It can tremendously simplify doctor's scanning procedure while improving patient's comfort.



Main features



- E-intelligence bed sheet replacing & automatic -correction system
- Ultrasonic coupling agent heater
- Vertical adjustment
- Heart position adjustment
- Emergency stop button
- The right-side bed flip function can support the patient in the side-lying position
- Trendelenburg/Reverse trendelenburg (optional)
- Horizontal adjustment (optional)



Product Functions

Heart position adjustment

The bed surface on the heart side can be flipped down to expand the inspection area during heart examination.

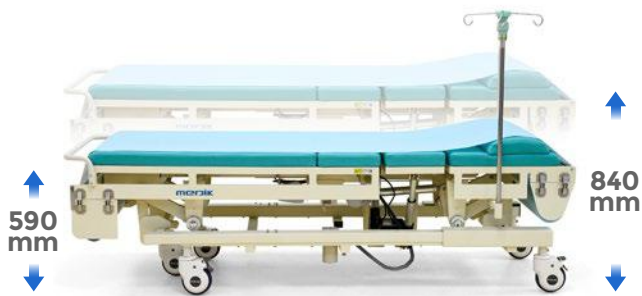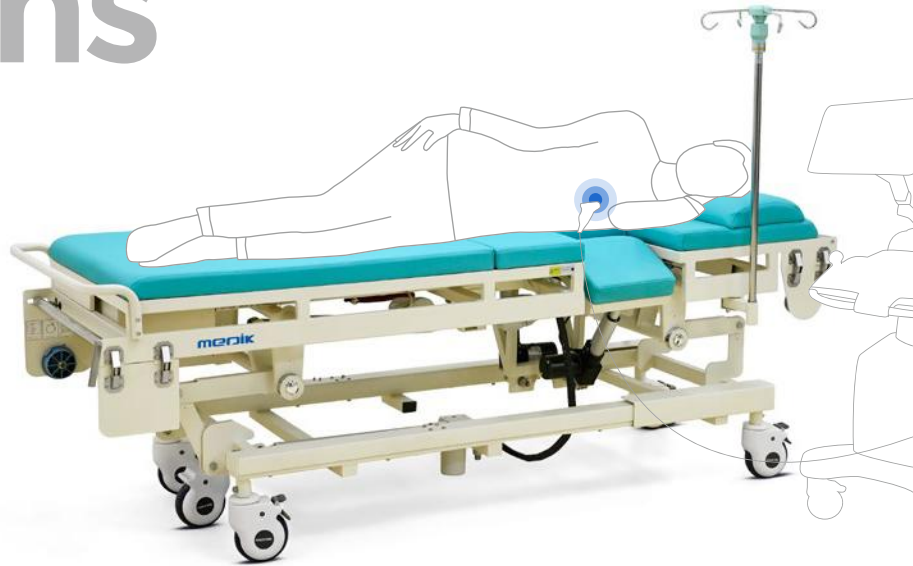

Height adjustment

590mm~840mm



Trendelenburg/ Reverse trendelenburg

0~12°



Horizontal sliding function

The bed surface has a horizontal sliding function, which is convenient for doctors to check



Back section adjustment

Convenient for patients to do lateral decubitus examination

Product Details



IV pole



Control button

There are control buttons at the end of the bed for the doctor's convenience.



Sheet roll holder

With automatic correction function and automatic paper changing function of three lengths.



Ultrasonic coupling agent heater

With heating function, can slide, easy to use.



Emergency stop button

Emergency stop button under foot section which can stop all movements urgently.



Hand remote

Infrared wireless remote control, easy to use



Castors

Five units of dia 5 inches castors with individual brakes.

Technical Parameter

Overall size	1900 * 650 * 590-840mm
Bed height adjustment	590mm~840mm
Trendelenburg/Reverse trendelenburg	0 ~12°
Horizontal movement	0~420mm
Load capacity	175Kg
Whole bed weight	133Kg

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US