

YA-H5-1

Manual Bed



*Better care
Better comfort!*

YA-H5-1

MANUAL BED

DESCRIPTION


The YA-H5-1 is a hydraulic bed, with effortless adjustment for height and pneumatic adjusting Trendelenburg. The heavy duty structure and easy adjustment is an ideal option for the area with power shortage.

FEATURES

- 1.Bed Frame electrophoresis and powder coating
- 2.Mesh bed platform, easy to ventilate
- 3.Stainless steel with MDF headboard and footboard, removable
- 4.back section adjusted by crank
- 5.Knee section is adjusted by crank
- 6.Height adjusted by hydraulic pump
- 7.Trendelenburg and Reverse Trendelenburg adjust by gas spring
- 8.Castors 125mm diameter ,2 guide and 2 with brake

STANDARD SPECIFICATION

Overall Size	2100*950*500-850mm
Lying surface	1950*900mm
Back section adjustable	0° to 60°
Foot section adjustable	0° to 40°
Height adjustable	500-850mm
Trend&reverse trend	0°to 10°
Safe working load	200kg

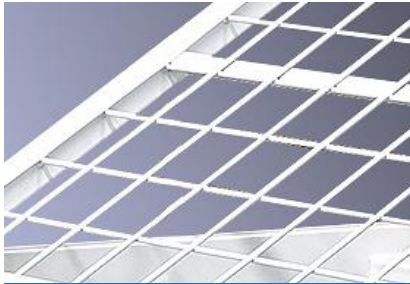

Mesh platform

Mesh platform provides great ventilation features to the patient especially when a long time in bed is needed.



Foot pedal

Height adjustment can be controlled by foot pedal.



Crank

Foldable metal cranks for back and leg adjustments.



Trend&reverse trend

The trend and reverse trend function is achieved by gas springs.



Hydraulic pump

Height adjustment is via hydraulic pump, durable and with large stroke.



Castor

Dia 125mm castors, 2pcs with brake and 2pcs with guide function.



Side rail

With aluminum alloy side rail to prevent patients from falling.



Head/foot board

stainless steel head/foot board with laminated panel gives an attractive appearance and easy to push the bed.



IV pole hole

There are 4 infusion holes around the bed, which is convenient for placing the IV pole.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US