

YA-M3-2

Manual Hospital Bed



*Better care
Better comfort!*

YA-M3-2

/ MANUAL HOSPITAL BED

The YA-M3-2 is a manual hospital bed with ergonomic and compact design. It is a unique manual bed with CPR function, this makes it suitable for ICU room with very limited budget.

Main Features

- Bed Frame electrophoresis and powder coating
- Slatted bed platform, strong and easy to ventilate
- Foldable aluminum side rail with 6 bar
- Stainless steel with MDF Headboard and footboard, removable
- 3 foldable cranks
- Castors 125mm diameter with diagonal locking

Technical parameter

Model	YA-M3-2
Overall Size	2160*940*500-750mm
Lying surface	1995*870mm
Back section adjustable	0° to 60°
Leg section adjustable	0° to 40°
Height adjustable	500-750mm
Safe working load	195kg

Main details




Slatted Platform

The slatted platform is strong and durable with good ventilation.



Crank

The foldable crank can prevent the handle from protruding outside to avoid collision.

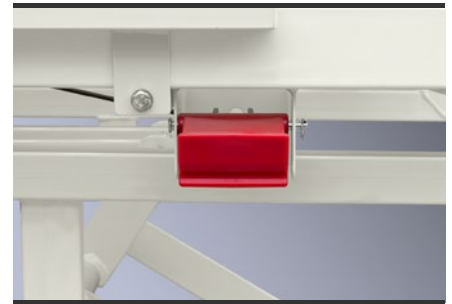

Manual CPR

Manual CPR handle at the head of the bed, which can quickly release the back section of the bed in an emergency.



Castor

Dia 125mm castors, each with independent brake.



Head/foot board

stainless steel head/foot board with laminated middle panel gives an attractive appearance and easy to push the bed.



Side rail

Side rails located on both sides of the bed guarantee the safety of the patient, and it is foldable down.



IV pole hole

There are 4 infusion holes around the bed, which is convenient for placing the IV pole.

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.