


# Manual bed

## YA-M5-8

This 5 function hospital manual bed with CPR function, It provide the high quality of construction, simplicity of use and safety. this manual bed is your best value options for better caring in ward room.

### MAIN FEATURES

- ABS head and footboard with holes for moving the bed easily
- Bumper castors on 4 corners for protection
- Dia 125mm Castors with brakes
- Four split side rails, drop down under platform
- Gas spring for back section adjustment and CPR function

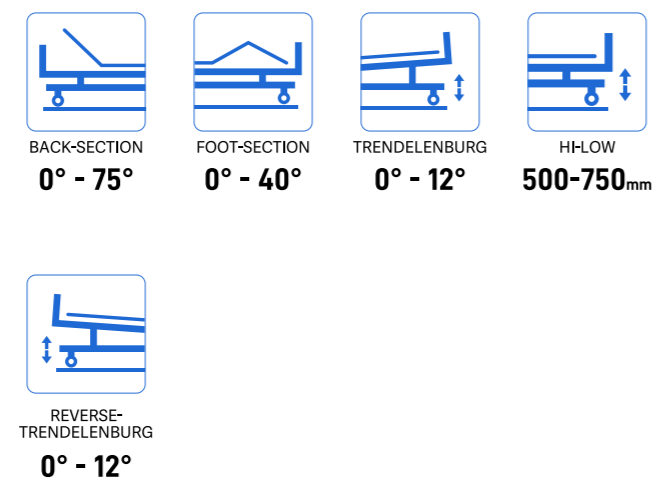

**CPR**  
Release the back-rest for emergency CPR.



**MANUAL CRANK**  
For manually operating back, knee, Hi-Low Trendelenburg & Reverse-Trendelenburg functions.



**CASTOR**  
Diameter of 125mm, Single castor with brake. (Option: castors with the central braking system)



**BUMP CASTOR**  
Reduce the impact of the collision, safe and slinky.



**CASE FOLDER**  
There is case card folder on the footboard to place patient information.



**HEADBOARD & FOOTBOARD LOCKING KEY**  
They can be easily fixed to the bed frame through this lock key and also easy to disassemble.



ITEM	PARAMETER
1	Overall size 2100*960*500-750mm
2	Backrest adjustable 0-75°±5°
3	Kneerest adjustable 0-40°±5°
4	HI-LOW 500-750mm
5	Trendelenburg 0-12°
6	Reverse-Trendelenburg 0-12°
7	Castors Four Dia 125mm castors
8	Loading capacity 220kg



No.3 Jinxiu Road, Jinfeng Industrial Zone, Jinfeng Town, Zhangjiagang, Jiangsu Province, China

+86 13962231292  
info@medik-medical.com

