



Pediatric bed

YA-PD5-2



Height adjustment



Section adjustment



Reverse Trendelenburg



Trendelenburg

MODEL:YA-PD5-2

PRODUCT INFORMATION

Pediatric bed

Description

Safety, attractive appearance, open and transparent environment are the most important features for a pediatric hospital bed.

The YA-PD5-2 with two column structure to minimize the entrapment risk also easy for hygiene, elevated ABS railing and optional cartoon design not only provide great safety for children but also make them feel peaceful on bed.

These child-friendly designs make it an ideal bed for PICU and pediatric care.



- The double-column structure can increase safety and facilitate cleaning
- Manual CPR and electronic CPR
- Auto-regress to remove the pressure of the abdomen when adjusting the backrest
- Auto-counter to provide the patient from sliding onto foot end
- Foldable side rails and removable head/foot boards made of ABS
- Central brake system
- Bumpers on 4 corners
- Urine bag holder
- Emergency button

MAIN FEATURES

Pediatric bed
MODEL:YA-PD5-2

option

- Nurse control
- Bed frame extension
- Weighing scale
- Night light under the bed
- Brake warning system
- 5th direction castor
- Battery backup
- Built-in control system of side rail


MODEL:YA-PD5-2

REMOTE CONTROL SYSTEM

Pediatric bed


1 Nurse control

Nurse control with multiple functions and locking button.

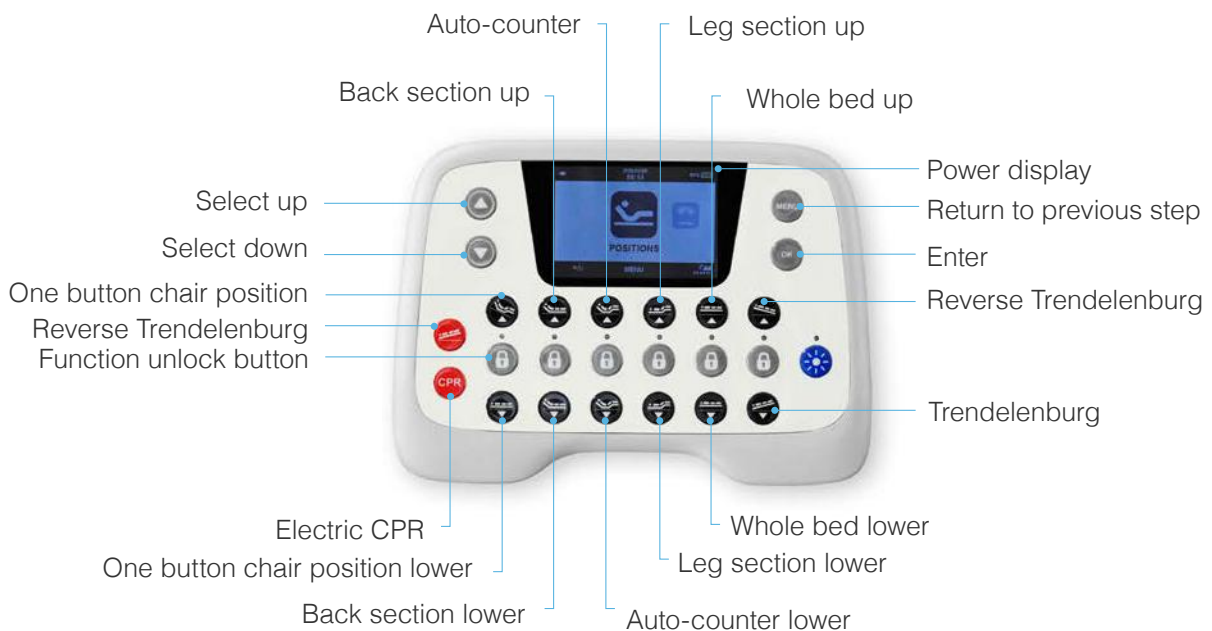

2 Hand remote control

The remote control has a lock function to prevent accidental touches.

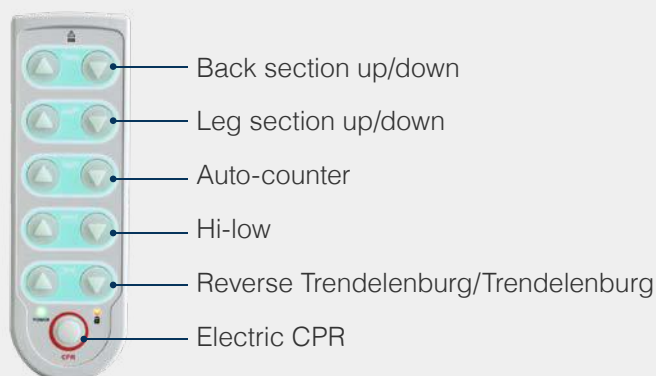

MODEL:YA-PD5-2

MAIN DETAILS

Pediatric bed

01



COLUMN STRUCTURE

The double-column structure can increase safety, prevent children from being pinched, and facilitate cleaning.

02



FOOT SECTION

Foot section can be lifted by gear.

03



IV POLE

Four provision holes for the IV pole around the bed.

04



SIDE RAIL

Four pieces of side rails with damper ,tuck-away type,adjusted through the buckle.

MODEL:YA-PD5-2

MAIN DETAILS

Pediatric bed

05



WEIGHING INDICATOR

There is a weighing display at the end of the bed, and it has a leaving bed alarm function.

06



CONTROL PANEL

There is a built-in control panel on the guardrail at the end of the bed.

07



ANGLE DISPLAY

08



BED EXTENSION

The bed extension is adjusted by a knob.

MODEL:YA-PD5-2

MAIN DETAILS

Pediatric bed

09



REMOVABLE HEAD AND FOOT PANEL WITH LOCKING DEVICE

10



NO BRAKE WARNING FUNCTION

When the bed stops but the brakes are not applied, the alarm will make sounds.

11



BUMPER CASTORS

Four bumper castors to protect the bed during transportation

12



MANUAL CPR FUNCTION

Manual CPR function on both sides of the bed for emergency situation.

MODEL:YA-PD5-2

MAIN DETAILS

Pediatric bed

13



EMERGENCY BUTTON WHICH CAN STOP ALL MOVEMENTS

14



URINE BAG HOOK

15



X-RAY BOX

X-ray box on the back section for easy access.

16



CASTORS

Four dia 150 mm castors , central brake pedal on both sides of the bed.

MODEL:YA-PD5-2

TECHNICAL PARAMETER

Pediatric bed

Model	YA-PD5-2
Overall size	1840*980*600-1000mm
Mattress size	1580*730*100mm
Bed extension	200mm
Knee section adjustment	0-30°
Back section adjustment	0-70°
Foot section adjustment	0-30°
Reverse Trendelenburg /Trendelenburg	0-12°
Loading capacity	240kg(530lbs)
Castor	Dia 150mm



Address : No.3 Jinxiu Road,Jinfeng Industrial Zone,Jinfeng
Town,Zhangjiagang,Jiangsu Province,China
Tel : 139 6223 1292
E-mail : info@medik-medical.com