

YA-PD5-4

Electric Pediatric Bed



*Better care
Better comfort!*

Electric Pediatric Bed

Model:YA-PD5-4




As a new generation of MEDIK's 5-function electric pediatric bed, YA-PD5-4 combines safety, attractive appearance, functionality and convenient operation. This model is widely accepted by the hospital and the children.

Safety

Children's safety are the most priority for the pediatric bed. The unconscious control and incapability should be taken into consideration and these are complete different comparing with an adult bed.



YA-PD5-4 is adopting 360° surrounding protection for the patient, and the height from platform to edge of railing can reach 500mm to provide enough safety from falling down even standing on bed.

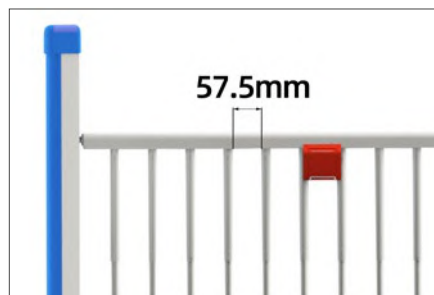

The gap between the railing is minimized to under 60mm to avoid entrapment of both heads and hands.



Double security railing locking mechanism located facing outside to deny the accessibility of children's mis-operation.

Main details



Attractive environment

Children are sensitive to the environment. To comfort the motion, YA-PD5-4i is designed with 360° visibility both from railing and the transparent head/foot board. Furthermore, the optional cartoon picture will also make the kids feeling relax.



Easy bed-exit

The unique telescopic railing provides not only safety but also easy bed-exit. The railing can be dropping down into a lower position to enable the children to get of the bed instead of being carried over from caregiver.



Functionality

As a medical bed which can be used in private ward or even ICU, it is featured with multi-functions. Motorized with four actuators, it has back, knee, height, trendelenburg and reverse trendelenburg bed positions.



Mobility

YA-PD5-4i is equipped with central brake, so one step pedal operation can unlock the bed or lock the bed firmly.

Main Features

- Full protection with elevated high railing for kids, 360° full enclosed protection
- Two ergonomics handrails make the release of the side rails extremely safe
- Transparent acrylic head and foot boards allow for good observation of the child's condition
- Central brake system with 150mm diameter dual-sided castors
- Fashion cartoon photos optional for head and footboard
- Telescopic siderail for easy bed exit

Technical parameters

Overall size	2000*850*480~780mm
Lying surface size	1820*720mm
Bed platform height	480~780mm
Trendelenburg/Reverse Trendelenburg	0-12°
Head and foot board length	630cm
Back section angle	0~75°
Foot section angle	0~40°
side rail height	500mm

ABOUT US

Zhangjiagang Medi Medical Equipment Co.,Ltd, founded in 1995, with more than 20 years experience in developing hospital furniture, is one of the largest manufacturer of hospital bed, patient stretcher, medical cart and other related hospital furniture in China or even Asia.

Through scientific management complying with ISO9001&ISO13485, by adopting advanced facility like Germany TRUMPF CNC, Japan OTC Robot Welding, Auto-power coating line with Swiss GEMA gun , Medi focuses on providing brand quality product with CE approval.

Headquartered in Zhangjiagang nearby Shanghai, with subsidiary oversea offices in Hungary, Dubai, Thailand, New Zealand, Kenya, Colombia, Medi acting professionally and swiftly for both pre-sales and after-sales.The products are widely accepted by global clients, such as UK RIE hospital, UNIVERZITETNI KLINICNI CENTER, Kokilaben Dhirubhai Ambani Hospital, Singapore National University Hospital etc.



YOU TUBE



MEDIK

Pls scan the QR code for more info

info@medik-medical.com | +86 139 6223 1292

Medik reserves the right to make changes without notice in design, specifications and models. The only warranty Medik makes is the express written warranty extended on the sale or rental of its products.

© 2022 Medik Services, Inc. ALL RIGHTS RESERVED. APR53703 rev 1 01-JAN-2022 ENG – US